

| | | | |
|--|----------------------|---------------|-------------------|
| | TECHNICAL DATA SHEET | 409701en | Date: 02/02/15 |
| | | Revision 7 | Page: 1/2 |

Codes: 409701, 409702, 409726, 409726.G, 409752, 409756 and 409756.B
150ml screw cap container (57 x 73 mm)

2014 catalogue pages: 154-155

| Code | description | sterile | case quantity | case weight (kg) | case volume (m3) | cases per pallet |
|----------|--------------------------------------------------|-----------|---------------|------------------|------------------|------------------|
| 409701 | Polypropylene labelled | no | 450 | 7.95 | 0.140 | 16 |
| 409702 | Polypropylene | no | 450 | 7.80 | 0.140 | 16 |
| 409726 | Polypropylene individually wrapped | Sterile A | 350 | 6.50 | 0.140 | 16 |
| 409726.G | Polypropylene individually wrapped | Radiation | 350 | 6.50 | 0.140 | 16 |
| 409752 | Polypropylene with PS spoon | no | 450 | 7.95 | 0.140 | 16 |
| 409756 | Polypropylene with PS spoon individually wrapped | no | 350 | 6.50 | 0.140 | 16 |
| 409756.B | Polypropylene with PS spoon individually wrapped | Radiation | 350 | 6.50 | 0.140 | 16 |

Directive 98/79/EC. "In vitro" diagnostic medical devices

| | No spoon | With spoon |
|------------------------|---------------------------|---------------------------|
| External diameter | 56.3 mm | 56.3 mm |
| Ext. diameter with cap | 62.1 mm | 62.1 mm |
| Total length (C) | 74.3 mm | 74.3 mm |
| Maximum volume (D) | 150 ml | 150 ml |
| Recommended volume (E) | 125 ml | 125 ml |
| Graduations | 20 / 40 / 60 / 80 / 100ml | 20 / 40 / 60 / 80 / 100ml |

Dimensions (±0.5 mm) Volumes (± 2 ml)



| Elaborated | Verified | Authorised |
|------------|----------|------------|
| | | |

| | | | |
|--|----------------------|---------------|-------------------|
| | TECHNICAL DATA SHEET | 409701en | Date: 02/02/15 |
| | | Revision 7 | Page: 2/2 |

SPECIFIC USE AND DESCRIPTION

Containers made of transparent polypropylene designed for collection of liquid or solid samples and transport to the laboratory for further analysis.

Each container features a discontinuous thread and a red ribbed cap made of high density polyethylene. High- security and leakproof seal. Caps ensure a positive hermetic seal from outside as well as from inside due to their inner ridge that perfectly fits the rim of the containers.


Containers bear a 3.3 x 5 cm frosted panel suitable for conventional markers. The panel has the following moulded-in inscription: Name, N^o, Date.

Code **409701** (labelled container): The label is designed to write the patient's name, the ward and its number, specimen type, time and date. Dimensions of the label: 55 x 37 mm.

Codes **409752**, **409756**, **409756.B** are suitable for the transport of faeces or other solid samples from the collection site up to the laboratory for analysis. The cap includes a laterally positioned removable spoon to facilitate solid sample collection. The removable spoon is 6.9 cm long.

Shelf life: 60 months for non sterile products. For Sterile A or Sterile R models: 48 months.

Storage recommendations: no special conditions, although sudden changes of temperature should be avoided.

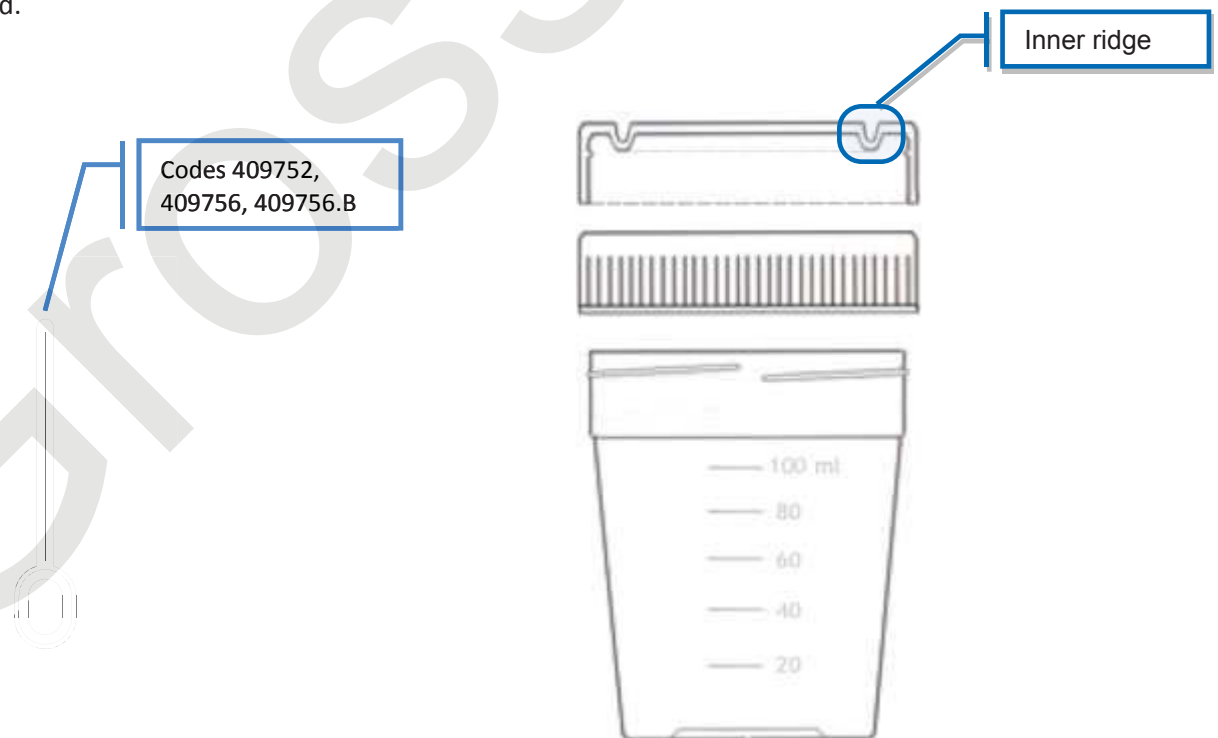

| Elaborated | Verified | Authorised |
|--------------|---------------|------------|
| <i>Gemma</i> | <i>Amalia</i> | <i>KEI</i> |